



BruxZir® Anterior Indications

- Individual anterior and posterior crowns
- Anterior veneers, inlays and onlays
- 3-unit anterior bridges

BruxZir Anterior Preparation Guidelines

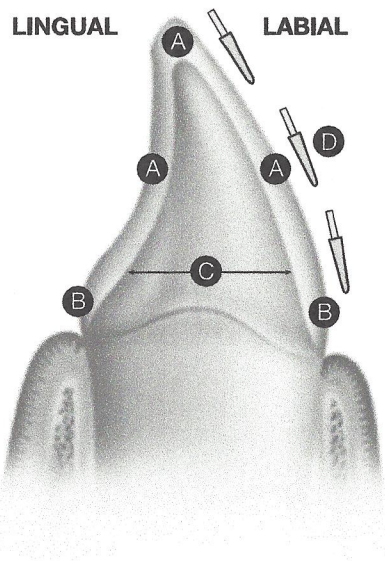

- A. 1.25 mm ideal reduction (0.8 mm minimum)
- B. Chamfer or shoulder margins preferred
- C. Labial and lingual walls must be convergent
- D. Preparation should be cut in three planes

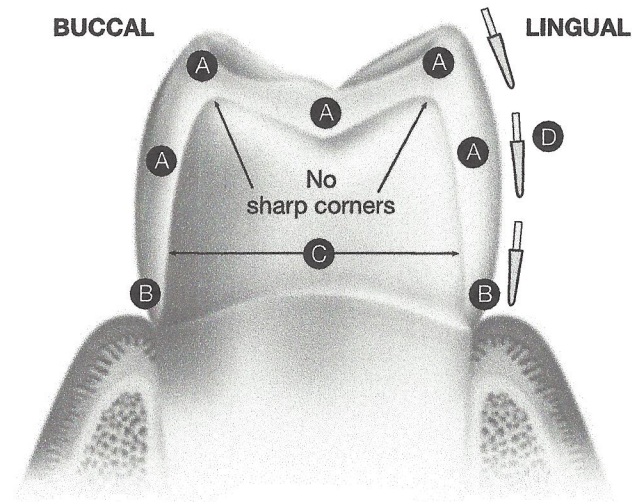

- A. 1.25 mm ideal reduction (0.8 mm minimum)
- B. Chamfer or shoulder margins preferred
- C. Buccal and lingual walls must be convergent
- D. Preparation should be cut in three planes